

BELL-HORN Premium Plus

since 1842

Code L0631 Applied

Ergo Premium Plus Instructions

Fasten the handle near the pulley unit. The handle is always on the left.

Stand holding the soft side panels apart, centering the back panel on the spine.

Wrap the left side to the front positioning the front panel. Wrap the right side on top of the left panel to close and attach in front. Adjust front or back panels as needed.

Pull the handle out and away from the power unit to desired compression and support. Place the handle on the fabric fastener band securely.

Tips:


Be sure to start with the handle near the pulley unit.

Extend or pull out the brace before application.

Adjust tension for sitting or standing, experiment to meet your needs.


Hand wash in warm water, rinse well, towel and air dry.

Ergo-panel, dynamics

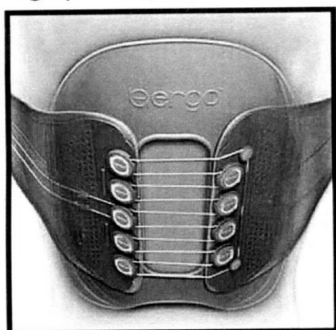


All Ergo-dynamic panels conform to each patient's individual body structure for improved comfort and support.

Ergo-panel height-adjustment

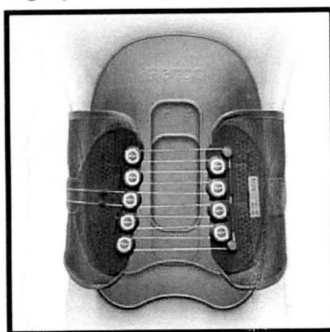


Standard Ergo-dynamic panels have 2 height settings adjustment--for additional patient comfort and support.

Ergo-panel, anterior attachment

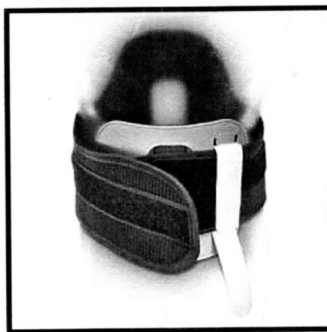

Ergo-dynamic anterior panels are easily attached with integrated velcro strap.

Ergo-dynamic panel easy and secure attachment (and removal)

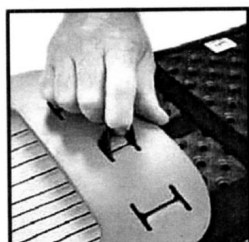


1. lift up attachment strap




2. place standard rivet head in slot



3. snap into place in slot

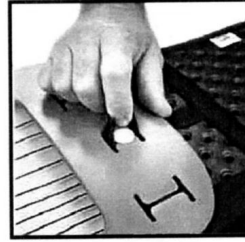

(placed rivet head)

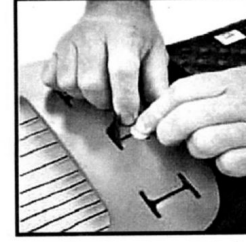

Slide rivet head to end of slot to remove